




DR and CR Protectors

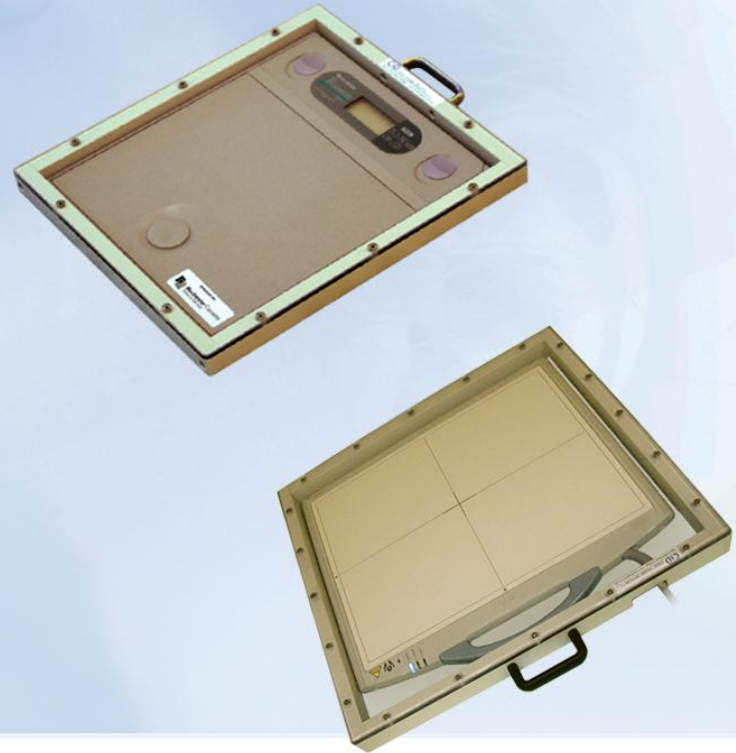

Keep Your DR Panels Safe and Avoid Breaking Your CR Cassettes!

You just purchased a top-of-the-line Radiography system and now you need to take weight-bearing X-rays of feet (AP view). The manufacturer says a patient can stand on the panel – but if it breaks, you're the one who has to pay to have it fixed or replaced. Ouch!

The Panel Protector provides a safe, easy way to take weight-bearing X-rays of feet (AP view). To use, simply place the unit over the imaging panel that is positioned on the floor. Patients step onto the Panel Protector with as much weight as needed to get the desired image.

The Panel Protector works with all CR systems and most DR systems, including DR systems offered by Canon, GE Healthcare, Philips, Quantum, Siemens, Shimadzu, Toshiba, Trixell, and Varian.

- *Panel Protectors come in all Standard CR cassette and DR Panel sizes*
- *Weight capacities:
For CR Cassettes – 275 & 450 lb (125 & 204 kg)
For DR Panels – 500 lb (227 kg)*
- *X-ray lucent faceplate allows X-rays to pass through the Panel Protector with no significant absorption or scattering*
- *Keeps patient's foot as close to the target image plate as possible, ensuring high quality images with no significant magnification*
- *Non-slip rubber grips keep the Panel Protector securely in place on a hard, smooth floor*
- *Easy to keep clean, the Panel Protector is made from germ-resistant, non-porous polyethylene, polycarbonate, and stainless steel*

Clear Image Devices ^{LLC}

3930 Michael Rd. North, Ann Arbor, MI 48103

Email: sales@ClearImageDevices.com

Tel.: 734.645.6459 Fax: 734.661.0220



www.ClearImageDevices.com

Copyright © 2012 Clear Image Devices LLC. All rights reserved.